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## High-pressure Preparative Chromatography Pump Installation manual

成都格萊精密儀器有限公司(Chengdu GeLai Technology Co., Ltd)  
太陽系分析基礎開發工作室(The Analytical Based Development Center)

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The GLP3500-500ml high pressure chromatography pump apply the latest ultra-high speed 32-bit RISC computing central processor + 32 bit DSP stepping driver processor as the core controller.

It is based on the frequency generation system, advanced motor drive system and friendly user interface that make the operation and control of the device very convenient.

Advanced pump head and check valve design ideas, a unique blend of the world's most advanced experience, allowing users in a variety of conditions , can guarantee the accuracy and repeatability specifications of infusion. The GLP3500-500ml dual piston high-pressure chromatography pump adopts dual pistons for one independent solvent delivery control.

The pump can be configured in isocratic or linear gradient mode.

While working in high-pressure gradient mode, the host computer or Master controller can control multiple pumps of the same type, and the gradient curve and gradient mode can be monitored through host computer for authentication.

Multiple communication modes, as the standard RS232, RS485, and Ethernet are selectable for pump system control.

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GLP 3500-500ml high pressure chromatogram pump(s) can be configured according to the needs of users to the following modes:

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High-pressure gradient Preparative	一個 2 元 至 4 元高壓泵系統. More than one GLP3500-500ml high-pressure chromatography pumps can be configured from binary(2 pumps) to quaternary (4 pumps) high-pressure gradient systems.
------------------------------------	---

### 1.3 主要特點 : Main features

泵頭 Pump head	根據需要，配備 500ml 泵頭 According demand to equip with a 500ml pump head
壓力自動校準 Automatic pressure calibration	使用者可以對顯示的壓力進行自動校準 Users can automatic calibrate pressure against displayed pressure
流量校準 Flow rate calibration	通過預設的 4 位元數位的改變，精確的校準該設備的輸液精度。 Accurately calibration of the flowrate precision of pump devices by the default four bits digital modulations.
人機界面 HMI	選用 192x64mm 的液晶顯示器，使得使用者編寫程式更直觀，更方便。 Adopt 192x64mm LCD monitor, so that the user may write programs more intuitively and more conveniently.
產品資訊提示 Product Information Tips	通過產品資訊的流覽，可以查詢儀器的出廠編號、安裝日期、泵的執行時間等資料，使用者可根據這些資料決定何時對該設備進行維護等。 By browsing through the product information, you can query the instrument serial number, date of installation, the pump execution time and other information, the user can decide when to carry out maintenance of the equipment based on these data.

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### The instrument principle and composition

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Introduction of the GLP3500-500ml High pressure chromatogram pump working principle.

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FIG.2 Block diagram of the principle of GLP3500-500ml

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The GLP3500-500ml high pressure chromatogram pump consists of the following five components:

- (1) 電源部件 power supply unit,
- (2) CPU 控制器單元 CPU controller unit,
- (3) 壓力測量控制單元 pressure measurement and control unit,
- (4) 馬達驅動器及馬達 motor drive and motor,
- (5) 泵主體元件部分 pump main element portion.

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Composed of control and drive power supplier, which provides powers of the 24V control and 75V drive that supplies to the pressure measurement and control unit, and the frequency generating unit for motor control.

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It is the core of the GLP3500-500ml high-pressure chromatography pump, which allows the user to set the required flow rates, to select the pump work mode, the edit time programs, etc.

### (3) 壓力測量控制單元 pressure measurement and control unit

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Measuring the size of the fluid pressure in the fluid path, the user can select a different values according to the column type and mobile phase; through the Max / Min pressure setting the high-pressure preparative pumps can be set to the optimum working condition, to avoid excess pressure and leakage.

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The power driver of the high pressure pumping system, which ensures liquid to flow through the column and flow cell steadily to have smooth preparation separations.

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Including components as the pump head and check valve elements, positioning screw mechanism etc.

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Check the environment requirements where the GLP3500-500ml is located , and should be met to:

供電電壓 Supply voltage	200~230VAC , 接地電阻 < 2Ω 200 ~ 230VAC, ground resistance <2Ω
環境溫度 Ambient temperature	20°C
環境濕度 Humidity	20%-80%
安裝場地 Installation site	安裝場所應遠離強震動源及強磁場干擾。 installation site should be away from strong vibration source and strong magnetic field

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According to user needs of different preparation demands to install 100ml, 250ml, 500ml, 1L or 2L pump head.

Pump head installation method:

Remove the two fastening screws of pump head, and remove the liquid line pipe that connect pump head and pressure sensor , install the required pump head to the pump body in the reverse order.

Change the appropriate settings in the pump settings display.

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After completing steps 3.2.1 and 3.2.2, connect the power cord into the instrument back panel and power on.

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Block diagram of the back panel of the GLP3500-500ml high pressure chromatogram pump

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Be sure to plug in the instrument power supply cord to three outlet pins power socket for correct source voltage, and good grounding.

Wrong power supply connections may cause serious personal injury.

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When first time to use the pump, run with small flow rate and to see if the fluid (mobile phase) can be pumped out. If it is not pumped out within a few minutes of continuous operation, stop pumping immediately and check the check valve and connecting tubes.

The pump head will be damaged when long-running with fluid (mobile phase) not sucked.

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All functions of the GLP 3500-500ml high pressure chromatography pump may be performed by 20 keys on the front panel.

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Figure 4 GLP3500-500ml high pressure chromatography pump schematic diagram

按鍵 Key	基本功能 Basic functions
-----------	----------------------

	 STOP	<p><b>返回鍵或停止鍵 Return key or stop key</b></p> <p>1、當液晶顯示器處於監控狀態時，該鍵為停止鍵。 When the LCD is in Monitoring state, the key used for Stop button,</p> <p>2、當在其它畫面時，該鍵為返回健，一旦按此鍵，將返回到監控畫面。 At other context, the key is pressed for Returning to the Monitoring state, When pressed, will return to the monitor screen.</p>
		<p><b>左移鍵或選擇鍵 Left or SELECT key</b></p> <p>1、當設置資料時，該鍵為左移鍵，左移被修改的資料位元，即閃爍顯示數位左移一位元。 When setting the data , the key is the Left Shift key to move the modified data byte left, namely move the flashing digit to left one position.</p> <p>2、產品維護控制時，此鍵為左泵頭控制選擇。 When used for product maintenance control, the key is for selecting the left pump head control .</p> <p>3、在正常監控畫面且泵處於未啟動狀態時，該鍵為一鍵清洗鍵。 In the normal Monitoring state and the pump is not activated, the key is used for One-Key Cleaning.</p>
		<p><b>右移鍵或選擇鍵 Right or SELECT key</b></p> <p>1、當設置資料時，該鍵為右移鍵，右移被修改的資料位元，即閃爍顯示數位右移一位元。 When setting the data , the key is the Right Shift key to move the modified data byte right, namely move the flashing digit to right one position.</p> <p>2、產品維護控制時，此鍵為右泵頭控制選擇。 When used for product maintenance control, the key is for selecting the right pump head control .</p> <p>3、在一鍵清洗鍵啟動後，該鍵為停止鍵。 It stop the running of One-Key Cleaning</p>
		<p><b>畫面鍵或遞增鍵 Previous Screen key or Increment key</b></p> <p>1、當按下該鍵時畫面翻轉到上頁。 Back to previous screen level.</p> <p>2、如果在資料設定狀態，按此鍵被修改位元的數字位元加 1，遞增範圍： 0—&gt;9—&gt;0 When setting the data, it is pressed to add 1 to the current value. Increment range: 0-&gt; 9-&gt; 0.</p>

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
	<p>When the LCD is in monitoring state, this key is Main menu key.</p> <p>A key / option in the main menu for returning to the top level, or Back to Previous button.</p> <p>When in parameter setting mode, the key is the CLEAR key. When pressed, the set value will be cleared.</p>
--	---

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Figure 5 GLP3500-500ml high pressure chromatography pump quick operation flowchart

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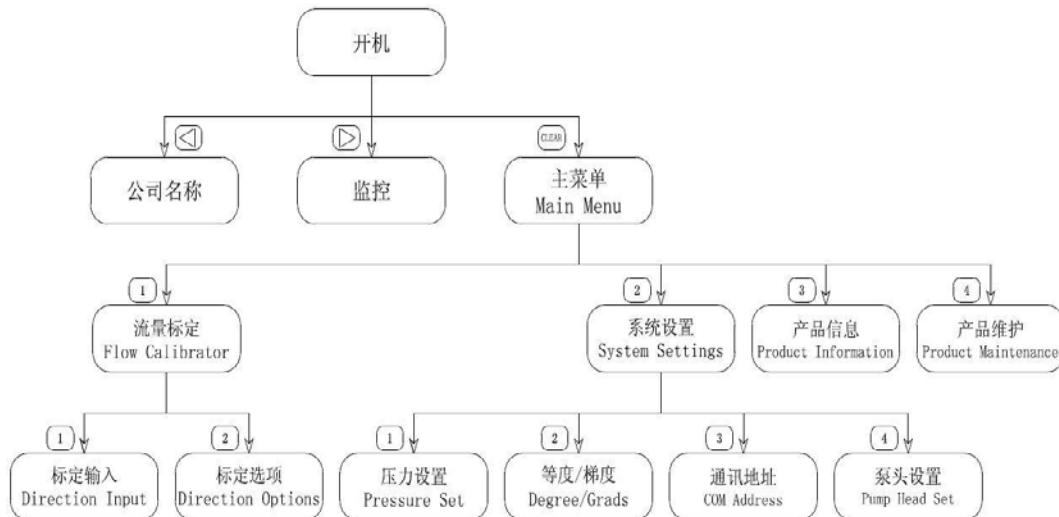Figure 6 GLP3500-500ml high pressure chromatography pump Begin Main-interface

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Figure 7 GLP3500-500ml high pressure chromatography Begin Sub-interface

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In the Main-interface, and press " ENTER/RUN " key to start the pump, press the "STOP" button to stop the pump. Press "HOLD" key to set the current acceleration time (Accel/min), working flow rate(mL / min).

Press the "Up" key to switch the screen as Begin Sub-interface, press the "Down" key to switch the Begin Main-interface screen, press the "CL" key to switch the main menu interface.

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In the main interface, press ENTER / RUN to start operation of the pump, press "STOP" button to stop the pump.

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Before starting the pump, you must ensure that the presence of liquid in the pump head.

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Before starting , run the system setup menu (System Settings) for settings, such as maximum, minimum pressure value (MAX, MIN), isocratic or gradient (Degree / Grads)

work mode selection, communication addresses (COM Settings) settings , pump head selection.

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In the main interface, press "HOLD" key to set the current acceleration time (Accel/min), working flow rate (mL / min).

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
It must to set a parameter value that is not too small or too large in their respective range, in order to Enter/Run to next parameter.

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
In the main interface, press "HOLD" key to switch the Main Menu interface as shown below:

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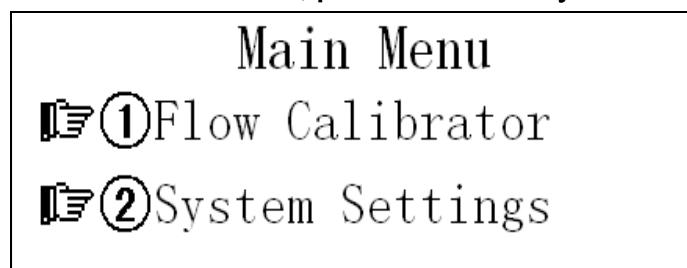
Option "①" digit key is for flow rate calibration (Flow Calibrator), press "①" digit button to enter the flow rate calibration interface. As shown below:

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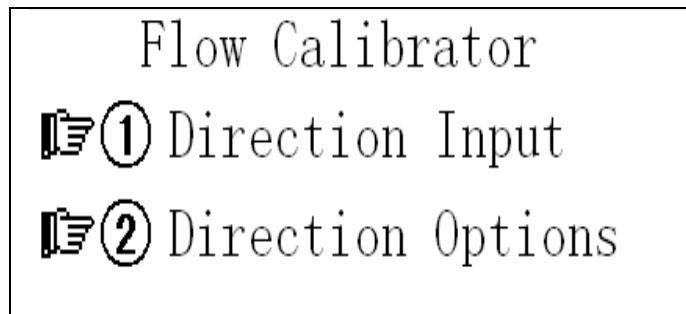In the between flow rate calibration sub-menu interface, press the "①" digit key to toggle to the calibration input (Direction Input) interface as shown below:

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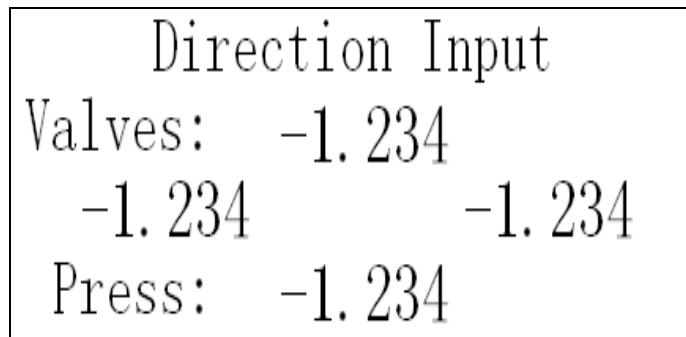Check valve coefficient (Valves), the pressure coefficient (Press)

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
In between flow calibration sub-menu interface, press "②" digit key to toggle to the calibration option (Direction Options)

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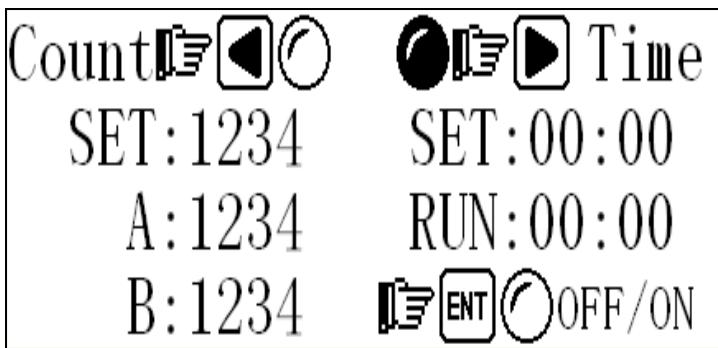The calibration mode selection: press the arrow keys "" for the counting mode; press the arrow keys "" for the timing mode; When you select the corresponding indicator should be highlighted.

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Calibration On/Off: Press "ENTER/Run" to toggle ON/OFF, with highlighted indicator corresponding to ON status.

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After the completion of the above setting, press "STOP" key to quickly return to the main control interface, press "ENTER" key to activate the pump to test flow rate.

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The calibration data has been completed in factory before the packging for shipping of the equipment , and is not required to modify under normal circumstances. As precision deviation caused by the wearing of equipment operation etc have been determined, it can be calibrated for modification.

If calibration parameter values updating is necessary, carry out in accordance with the following program:

**(例) 測常壓 (單向閥)****(example)Measuring normal pressure (check valve)****1、準備好量筒**

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**NOTE: The pump has a self-test function for piston position recovery at every each running start, and the start to work after the completion of the self-test.**

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After the completion of the self-test, immediately put the fluid outlet tube into the measuring cylinder, and the pump will start running.

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When the pumps start to run, the A pump head and B pump head begin to count, when the count value (time) equal to the set value, the pump automatically stops.

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Read the actual value of the liquid collected in Measuring Cylinder(180mL)

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coefficient = new calibration value,  $180/(50+50)-1.825=-0.025$

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Enter the calculated coefficients, and change to count by timing, then enter a time value (00:10) (10 minutes) to verify the accuracy of the new factor. Repeat steps 4-7, if the Actual value is correct, turn off the calibration switch and you can start to work. If there is an error, repeat steps 1-8 until the flow rate value satisfied.

Before to change flow rate calibration parameters (value), ensure that this value has been

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It will lead to very cumbersome recalibration, if the number is forgotten.

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
NOTE: The parameter settings of following menu is same as the previous Flow rate calibration

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
In the main interface, press "CLEAR" button to switch to the main menu (Main Menu) interface, as shown below:

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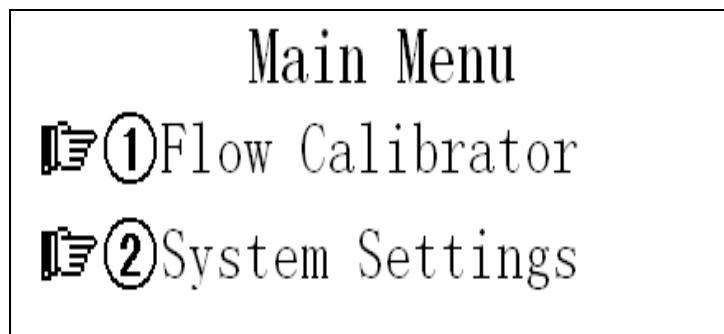The "②" digital key is for system settings (System Settings), as shown below:

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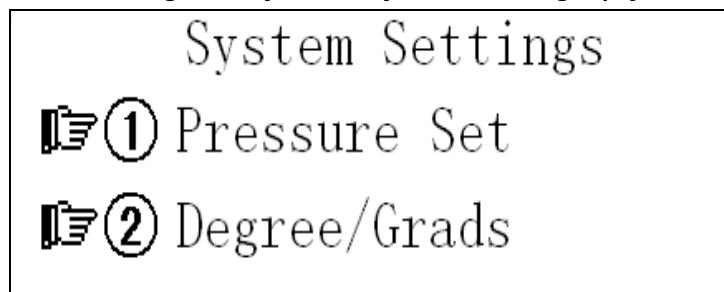In the System Settings (System Settings) interface, the options "①" digit key is for pressure setting (Pressure Set), as shown below:

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when the system pressure exceeds the set value, GLP3500-500ml will automatically stop and alarm alert tone.

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when the system pressure is maintained lower than the set value more than two minutes, GLP3500-500ml will automatically stop and alarm alert tone.

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In the System Settings (System Settings) interface, the options "②" digit key is for isocratic / gradient setting(Degree / Grads) , as shown below:

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Choose 0 for system to operate in isocratic mode.

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Choose 1 to operate in high-pressure gradient mode.

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In the System Settings (System Settings) interface, the options "③" digit key is for communication settings (COM Settings), as shown below:

COM Address	
Address: 000	
001=A	003=C
002=B	004=D

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Select 1,2,3,4 that correspond to pump A , B , C , D respectively.

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
The gradient curve reflects the mobile phase composition percentage variation of the pumps in actual operation.

SHOW GRADIENT can explain to the user to check the correctness of the gradient programming, whether the pump programs is running in consistent with the designer's ideas.

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
See the example gradient curve in Figure 8
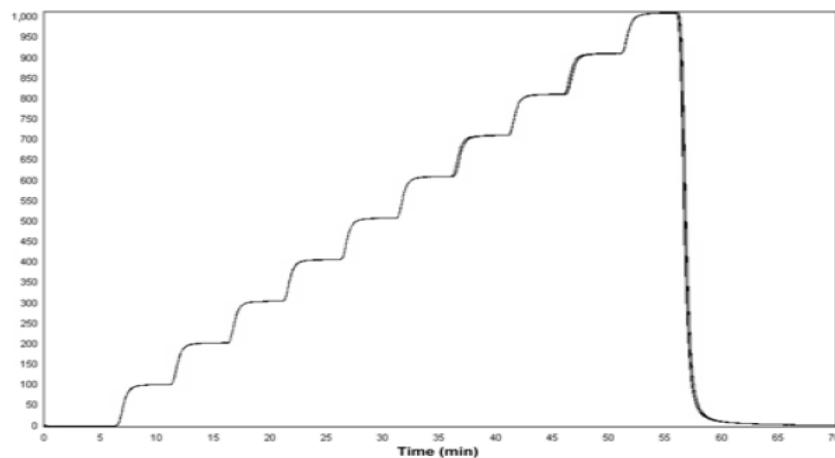
## 4.6 產品資訊 ( Product Information ) Product Information

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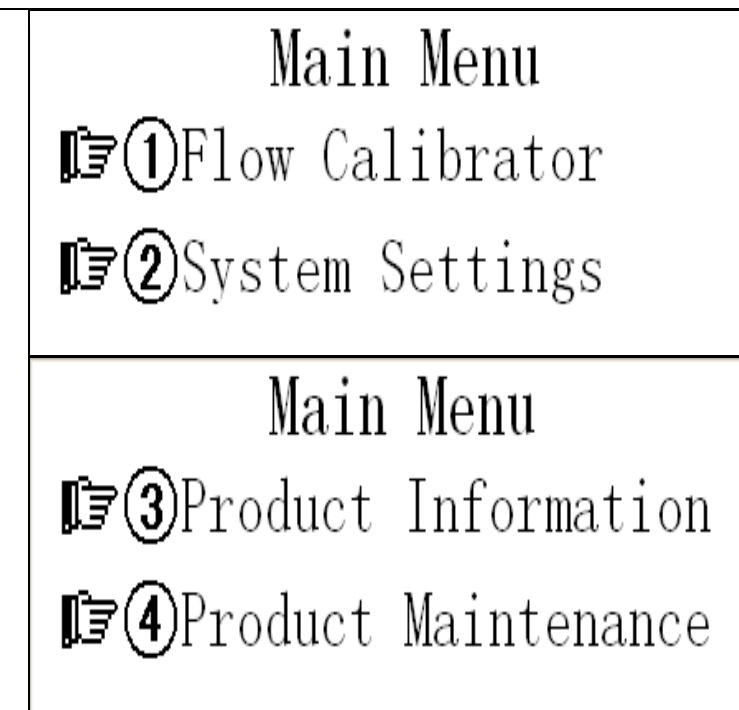Product Information
Run time: 1234567:00:00
Product NO.: GLP 123456789
Factory date: 0000-00-00

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
## 4.7 產品維護 ( Product Maintenance ) Product Maintenance

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
In the main interface, press "CLEAR" button to switch to the main menu (Main Menu) interface as shown below:

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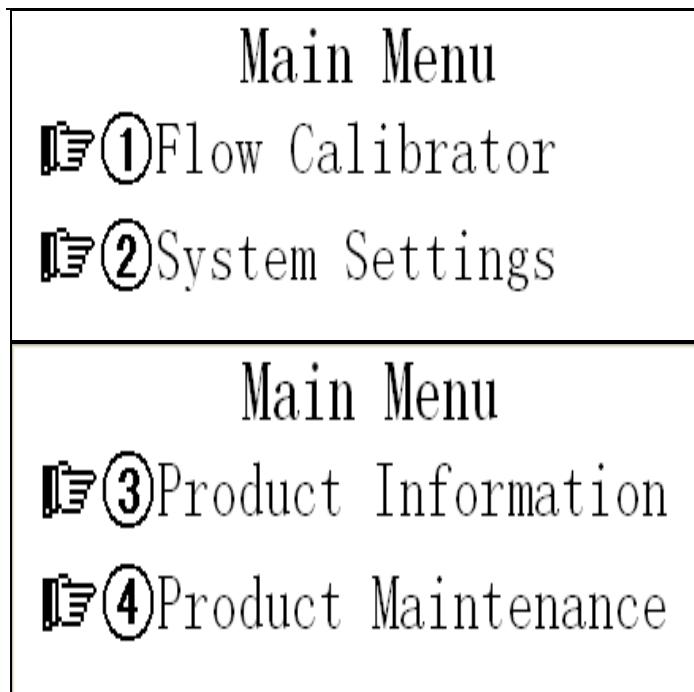Options "④" digit key is for Product Maintenance, as shown below:

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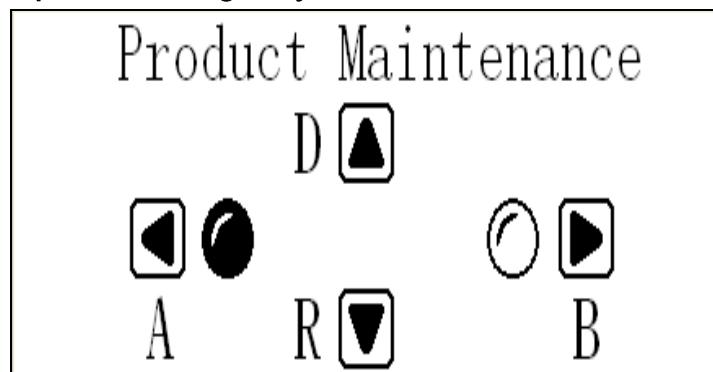
In the main interface, press "HOLD" key to switch to the Main Menu interface as shown below:

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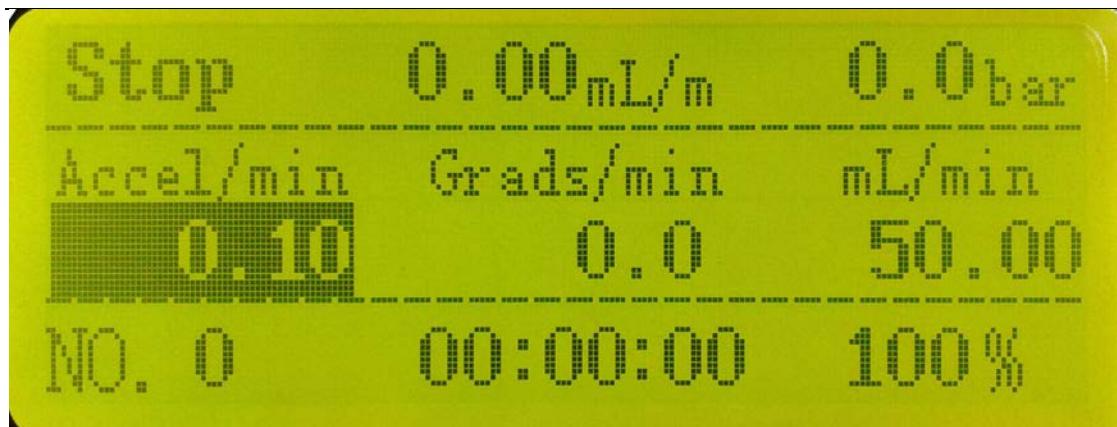Press the HOLD key to start modifying parameter settings, "Accel/min" is to modify the pump start acceleration time. The parameter currently being modified is highlighted, wherein the number of byte to be modified flashes. When the parameter is set to the desired values, press ENTER/RUN key to save the setting. Press HOLD again before ENTER/RUN key pressed If the modification of this parameter is not necessary to cancel the current modification, and skip to continue editing the setting of next parameters.

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High-pressure Preparative Chromatography Pump Installation & User Manual

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Press key to turn to down level menu, and press the operator #5 key to enter the flow rate setting interface

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Press the HOLD button to start modify the parameter setting value and then press ENTER/RUN to confirm the flow rate. Press STOP to return to the main control page, and press the key on the panel to start pump cleaning flow. The pump cleaning flow can only be stopped by pressing the key on the panel.

The STOP key can not be used to stop the pump cleaning flow.

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More than one GLP3500-500ml high-pressure chromatography pumps can be composed to binary to quaternary high-pressure gradient chromatography systems.(See figure)

The dynamic mixer (GLP3500-500ml-DM) is necessary to be configured here.

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Dedicated serial cable is necessary to connect the GLP3500-500ml high-pressure

chromatography pumps to the computer communications port.

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The serial ports connection diagram of High-pressure gradient system

## 六、簡易維修 Basic Maintenance

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Installation and replacement of the pump head

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1. Ensure the accomplishment of the pump head Installation, set up and adjustment.
2. Advance the pump head in the vertical direction
3. Rotate clockwise the two pump head installation screws 2-3 circles,

4. Tighten one pair of diagonal screws slowly.
5. Tighten the other pair of diagonal screws slowly, then for .
6. Confirm all the four screws have been tightened
7. The disassembly of the pump head is reversed to above process.

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If the pump head needs cleaning or corrupted, the above disassembly and installation steps should be applied.

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The check valve should be cleaned once it becomes dirty. If not cleaned , it may affects flow rate accuracy of the pump, and also can cause larger pressure fluctuations, and can result in failed infusion in serious cases.

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- (1)Remove the check valve from the check valve seat.(Note to memo the check valve direction to the check valve seat )
- (2)Carefully remove the check valve.
- (3)Place the check valve into the beaker using ultrasonic cleaning in acetone for 20 minutes.
- (4)Install the check valve back onto the check valve seat. (Note that the check valve direction must be consistent with that it is taken out from the check valve seat).

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Under the premise in accordance with the correct operation rules, we have 12 months warranty after the date of installation.

(except the damage to the sealing parts)

Material costs may be charged if damage is caused of human factors, such as damage improper operation of the instrument in all periods including the warranty period.

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GLP3500-500ml high pressure chromatography pump has been properly packaged at the factory. If it occurs damage caused by transportation , the customer should declare as soon as possible to the relevant departments to notify the manufacturer for repair or replacement.

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Infusion mode Screw drive parallel twin piston reciprocating pump, double piston pipette

最大工作壓力	Maximum operating pressure
GL3100	35MPa
GLP3250	20MPa
GLP3500	15MPa
GLP3000-1L	12Mpa
GLP3000-2L	10Mpa
流量範圍	Flow rate range
GL3100	0.1 ~ 100ml/min
GLP3250	1 ~ 250ml/min
GLP3500	1 ~ 500ml/min
GLP3000-1L	10 ~ 1000ml/min

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## Ordering & Service Information: 訂購和服務信息.

Welcome to our web site to see our currently updated product information and catalog.

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NO.137, JIAZIDONGDAO, DAYI COUNTY, CHENGDU, CHINA 611330

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電 話(Phone) : 028-88201806/88201866, 傳 真(Fax) : 028-88201866,

Web: <http://www.cdgelai.com>, Email: [cdsglp@163.com](mailto:cdsglp@163.com)

The Analytical Based Development Center (ABDC Work Shop)

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Zipcode:407, Add: 5th Floor, NO.641, Fu Shun Road, Shi-Tuen District , Taichung City, Taiwan, R.O.C.

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[http://chromnet.net/Product%20Atlas\\_Preparative%20Chromatography%20Systems.aspx](http://chromnet.net/Product%20Atlas_Preparative%20Chromatography%20Systems.aspx)

Email: [service@chromnet.net](mailto:service@chromnet.net)

LINE:service.abdc, SKYPE: skypeabdc, Twitter: ttabdc, WeChat: wcabdc

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ChengDu GeLai Technology Co.,Ltd.